

# AMRAD Medical



## DFMT ELITE

PREMIER RADIOGRAPHIC SYSTEM

# Performance, Versatility and Value

Designed for high volume throughput for hospitals, orthopedics and medical imaging facilities. The DFMT Elite can perform a broad range of exams while optimizing technologist productivity and patient workflow. Simplicity and functionality provide broad positioning capability for patient procedures, including cross table laterals. Feature-rich controls and ease of use enables the technologist to focus on patient care.



## Overhead Patient Handgrip

For lateral exams and weight bearing studies for 360° positioning with increments at every 5° for maximum flexibility

## Remote Collimator Light

Conveniently turn collimator light on from the wall stand

## Smart Button Technology

Rotating the angulation handle from table to wall stand, the buttons for vertical and horizontal movement intuitively swap functions for seamless operation

## Table Height Indicator

Ensures optimal working height

## Adaptive Optimal Display

As the angulation handle is rotated to the wall stand, the display changes orientation for easy viewing


## Simple Service

LED operational status indicators

## T.E.C.H.

Table Elevation Control Handle for table top lock release and vertical movements

# Features and Benefits

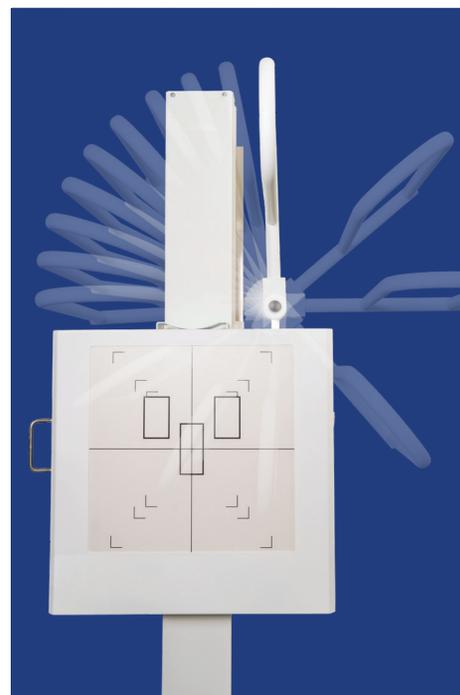

## Tube Stand

- Floor mounted rail for longitudinal positioning to table and wall stand
- Transverse and rotational positioning for broad range of patient exams
- Adaptive Optimal Display shows continuous SID to table and wall stand
- Smart button technology for ease of positioning

### Options:

Trunnion mount

## Table

- Elevating, four-way float, flat top for easy transfer and positioning of patients
- 650 pound patient capacity
- Recessed four pedal design provides easy operation and patient safety
- Collision protection and fail-safe locks for added patient safety
- Table Elevation Control Handle (T.E.C.H.) with lock release
- Indicator light for optimal working table height

### Options:

35.5 inch wide table top

## Wall Stand

- Center or side mounted for patient and/or configuration flexibility
- Ergonomic lock release handle for ease of positioning
- Side mounted and overhead hand grip for patient stability and positioning
- Remote collimator light improves productivity

### Options:

Dual sided lock release

Long format 17 x 51-inch

# Innovative Orthopedic Application



## 17 x 51-inch Long Format Wall Stand

- Single exposure produces 17 x 51-inch image
- Long leg, scoliosis, full spine exams
- Isolate top, middle, or bottom for single 17 x 17-inch image
- Dual lock release handles and dual remote collimator light features
- Available with new full rooms suites or as a retrofit to existing x-ray suites

## The DFMT Elite...where simplicity meets functionality

- DFMT Elite is fully featured and can be configured with a variety of generators and X-ray tubes to meet the clinical requirements of any imaging environment.
- Traditional floor mount design means no learning curve for the technologist.
- Industry leading Warranty, Customer Service and Technical Support.
- The DFMT Elite is compatible with all DR solutions.



[www.amradmedical.com](http://www.amradmedical.com)

SUMMIT INDUSTRIES, LLC  
7555 N. Caldwell Avenue, Niles, Illinois 60714  
(p) 773-353-4000 | (p) 1-800-729-9729



LITR038-1119