
# DEL MEDICAL OTC12 Overhead Tube Crane Systems

## Introducing the Entirely NEW

#### Introducing the Entirely NEW Del Medical OTC12 System

The Del Medical OTC12D system (pictured) is a high-throughput imaging suite offering complete versatility and reliability at an affordable price. The durable OTC12 overhead tube crane combined with our 800lb patient load elevating table and slender design wall stand can readily meet the demands of the busiest imaging environments with ease.



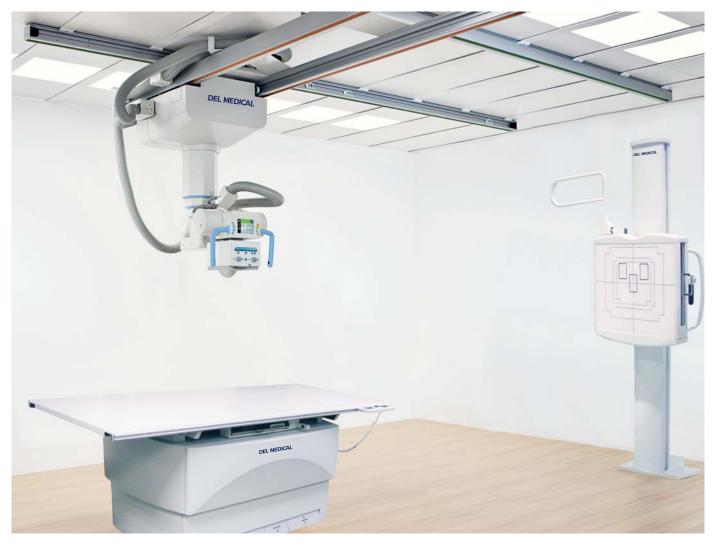



## DEL MEDICAL OTC12D Overhead Tube Crane System

#### OTC**12** Overhead Tube Crane

This overhead mounted support allows for the most flexible and affordable configuration choices in radiography. Sturdy and reliable, it is ideal for high-throughput environments such as hospitals, emergency rooms and high volume out patient facilities.

- Touch screen display for SID and angulation
- Telescoping arm allows for small focus-to-ceiling distance
- Ergonomic design with easy movements for fast, accurate positioning
- Electromagnetic brake for telescopic lift and transverse track
- Integrated lock function
- Available with manual or automatic collimator (pictured) with vertical tracking capability



#### EV**800**

Elevating, Four-Way Float Top Table

This rugged table is built to demanding hospital quality standards and features an 800lb (363kg) patient load capacity allowing complete confidence when imaging larger patients. The smooth rounded corners and integral collision sensors ensure patient and operator safety.


- 17"x17" bucky with 10:1 grid and manual or size-sensing cassette tray; an optional DR ready bucky is also available!
- Smooth and quiet table operation with rack and pinion motor drive
- Features recessed foot treadles and table-side hand control for easy control of table movements
- Smooth, stain resistant fiber resin table top with variable height from 22.25" (56.5cm) to 34.25" (87cm)

Available table options include multiple table top sizes, carbon fiber top, a variety of grids, patient hand grips, compression belt and lateral cassette holder.





Table side hand control provides additional source for all table movements

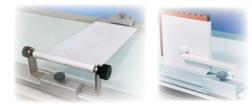

Optional compression band and lateral cassette holder



Table top lowers to 22.25"(56.5cm) allowing for safe and easy patient transfers


#### VS**300** Vertical Wall Stand

This robust wall stand features the same streamlined design with soft rounded edges that compliment the  $\mathsf{EV800}$  table.

- Slender column design with cable concealment system
- Electric "Fail Safe" locks and ergonomic release handle
- Patient handgrips, two lateral and one removeable overhead
- 17" x 17" bucky with 10:1 grid and manual or size-sensing cassette tray; an optional DR ready bucky is also available!
- Receptor lowers to floor to accommodate a full range of upright radiographic procedures
- · Counter balanced with smooth, vibration free movement

Available options include floor mount kit for free standing version and a variety of grids.

Extensive SID travel range from 72in (182.9cm) to 15.5in (39.4cm)

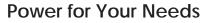

Choose from Del Medical's CM and CMDR Series of 50, 65 or 80kW APR generators. Compact in size and lightweight, these generators offer high frequency performance up to 400kHz and come complete with pedestal console and handswitch. Our CMDR generators have digital ready functionality on select DR systems so they are DR ready when you are ... now or in the future.

Combine with any of our 300, 400 or 600 kHU x-ray tubes to build a system with the versatility to fulfill the needs of demanding imaging suites.



A DELWORKS DR Upgrade Package is the easy and affordable way to take the next step to digital. Choose from a variety of fixed or portable, tethered or wireless, flat panel detector(s) complete with DELWORKS acquisition and processing workstation and software. Please contact us for complete details on our DR upgrade packages.





#### OTC**12S-A** Overhead Tube Crane System with Tilting Wall Bucky

The Del Medical OTC12S is our premier imaging suite that includes:

- OTC12-A with automatic collimator, vertical auto tracking and touch screen display for SID, angulation and generator controls
- Rad TC-A elevating table with 617 (280kg) patient load capacity
   and automatic size-sensing cassette tray
- WS 04 wall stand with tiltable bucky from -20° to +90° with 0° detent
- IT Series high frequency, anatomically programmed generator with Plus Package (heat calculator and AEC). Choose from 55, 65, or 80 kW

Contact us or your local Del Medical dealer for complete specifications!



EDICAL

## OTC12 Series Overhead Tube Crane





#### Specifications

- OTC12 Overhead Tube Mount
- Minimum focus ceiling distance:
- Vertical telescope travel range:
- Longitudinal travel range, -with standard rail:
   -with optional rail extensions:
   -detent positions:
- Transverse travel range, -with standard 9.8' (3m) rail: -with optional 13.1' (4m) rail: -detent positions:
- Transverse track length, ±39"(10mm): 119" min. (3030mm)
- Tube rotation range, horizontal axis: -detent positions:

**OTC12 Overhead Tube Mount Options** 

• 300kHU, 400kHU or 600kHU 0.6-1.2mm, X-ray tube

· Collimator (Ralco manual, OTC12M or Siemens

• Extended longitudinal rail length, up to 20' (6.37m)

Extended telescope and travel, increases vertical travel by 7.9" (20 cm)

- Tube rotation range, vertical axis: -detent positions:
- Maximum installation height:
- Touchscreen interface:

automatic, OTC12A)

Vertical tracking

32.7" (830mm) 70.8" (1800mm)

max. 136.2" (3460mm) 219.9" max. (5585mm) configurable during installation

86.6" min. (2200mm) 139.7" max. (3550mm) configurable during installation 119" min. (3030mm) 172" max. (4380mm) 120°, +120° -90°,0°, +90° +182°, -154° +180°, -90°, 0°, +90° ≤1.2 miles (2000m) above sea level SID and horizontal tube rotation. Generator control on OTC12-A-S suite only.

#### EV800 Elevating Table

Base length:	53" (135cm)
Base width:	25" (64cm)
Top length, standard:	86.5" (220cm)
-optional	78" (198cm)
Top width, standard:	31.9" (81cm)
-optional	35.9" (91cm)
Table weight:	563lb (256kg)

See Del Medical EV800 component spec sheet for complete specifications

#### VS300 Vertical Wall Stand

 Height:
 84" (213)


 Depth:
 13.4" (3-2)

 Width (with receptor)
 24.6" (6-2)

84" (213.4cm) 13.4" (34cm) 24.6" (62.5cm)

See Del Medical VS300 component spec sheet for complete specifications

Some flat panel DR detectors may have some affect on table and/or wall stand bucky specifications, depending on manufacturer. Please contact us for complete details on digital upgrades!

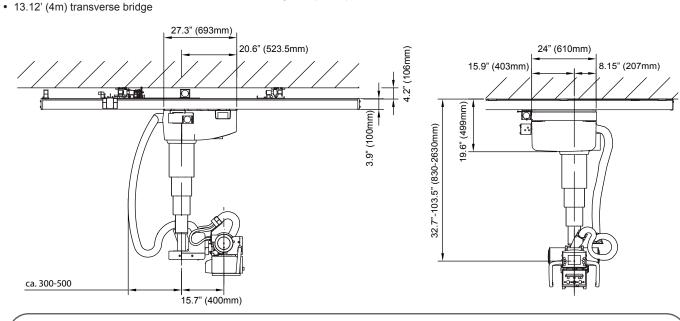

## DEL MEDICAL

Toll Free: (800) 800-6006Del MFax:(847) 288-7011SpectEmail: info@delmedical.comwww.delmedical.com241 Covington Drive, Bloomingdale, IL 60108

Del Medical is a UMG Company. Del Medical is an ISO 13485 & ISO9001 Certified Facility. Specifications are subject to change without notice.

Toll Free: (800) 261-9808 Fax: (914) 835-6111 Email: info@umgxray.com www.umgxray.com 28 Calvert Street, Harrison, NY 10528

05.13 #9000-DelOTC12